
INSTRUKCJA

OBSŁUGI

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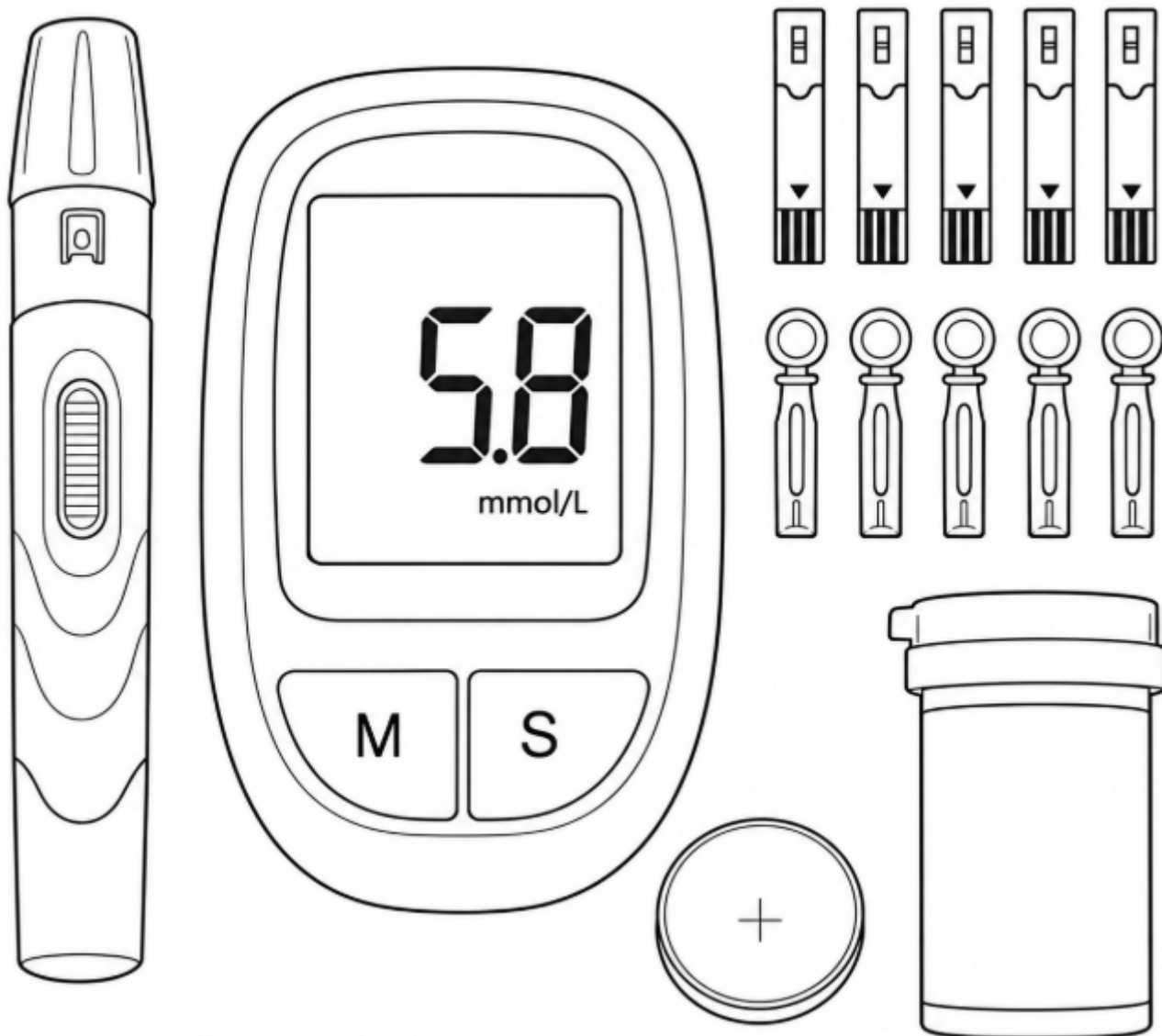Dear Sir or Madam, thank you for purchasing our product!

Before using the product, please read carefully with the instructions below for proper use of the product.

Please keep this manual for future reference and following its recommendations, because failure to follow it may pose a threat to life or health.

Introduction to the user manual

This user manual has been prepared to help you familiarize yourself with the proper and safe use of the device. It contains important instructions on operation, maintenance, and safety rules, which, if followed, will ensure the smooth operation of the device and extend its service life.

Before first use, please read the instructions carefully. The information contained herein will also help you properly store the device, especially if it will not be used for an extended period.

Following the recommendations below will allow you to enjoy trouble-free and safe operation of the product.

ATTENTION!

Using the device in a manner inconsistent with the information contained in the manual may result in the risk of bodily injury.

Application and description of the device

The meter is intended for self-monitoring of glucose levels in fresh capillary blood collected from the fingertip . According to the product data, measurements can also be performed on the forearm, provided that the glucose level is stable and in accordance with the device's operating instructions.

This device is intended for home use only and is not intended for medical or commercial use, and is not intended to replace medical therapy. The device is intended for indoor use in dry rooms. It is not intended for medical, therapeutic, or commercial use.

For safety and CE certification reasons, this product must not be rebuilt or modified in any way. Using the product for purposes other than those described may result in damage to the product. Improper use may also result in hazards such as short circuits, fire, electric shock, etc.

WARNING!

- Danger of improper use!
- Using the device in a manner other than its intended purpose and/or using the device in any other way may involve various risks.
- The device should only be used for its intended purpose.
- Please follow the procedures described in this instruction manual.
- Keep the product away from children.

TIP

- The device must be checked for completeness of delivery and any visible damage.
- In the event of an incomplete delivery or damage due to faulty packaging or transport, please contact the service hotline.

Technical data of the device

Symbol	27729
EAN	5907451382509
Material	plastic
Display	LCD
Power supply	Button battery
Measurement time	about 8 seconds
Result units	mmol/L or mg/dL
Memory	up to 128 results

PREPARING FOR FIRST USE

- Prepare your meter, test strip, lancing device and new lancet.
- Wash your hands with warm water and soap.
- Dry your hands thoroughly.
- Check the condition of the test strip.
- Do not use a dirty, wet or damaged strap.
- Remove the strip from the container immediately before taking a measurement.
- After removing the strip, close the container carefully.
- Place a new lancet in the lancing device.
- Do not use a previously used lancet.
- Make sure the meter and all accessories are clean and dry.

Startup/Use

- Please remove the packaging materials and dispose of them properly and clean the product and accessories.
- Insert the batteries as marked, paying attention to the correct polarity (+ and -).
- Close the compartment and make sure the battery is seated correctly.
- Wash and dry your hands thoroughly.
- Remove one test strip from the container.

-
- Insert the strip into the meter slot in the correct direction.
 - Once the strap is placed, the device will turn on automatically .
 - Wait until the screen displays a message indicating that you are ready to collect a blood sample.
 - Prepare the lancing device with a new lancet.
 - Place the lancing device on your fingertip.
 - Make a puncture.
 - Obtain a small drop of fresh capillary blood.
 - Bring the front edge of the test strip close to the drop of blood.
 - Allow the strip to sample according to its design.
 - Do not smear blood on the strip and do not press your finger directly onto the surface of the strip.
 - Once the sample has been collected correctly, the device will start measuring.
 - Wait about 8 seconds .
 - Read the result from the LCD display.
 - Check what unit the result was presented in: mmol/L or mg/dL .
 - The result will be saved in the device's memory.
 - Once the measurement is complete, remove the test strip.
 - Do not reuse a used strap.
 - Remove the used lancet from the lancing device carefully.
 - Dispose of used strips and lancets safely.
 - After completing the test, wash your hands if they have come into contact with blood.

RESULTS MEMORY

- The meter can store up to 128 results .
- Each record can contain the measurement result and the corresponding date and time.
- The device buttons are used to view the saved results.
- The memory function allows you to compare previous measurements.
- When using the memory, pay attention to the date, time and marking of the given result.
- The device's memory is intended to assist in self-monitoring and does not replace medical documentation or interpretation by a specialist.

AVERAGE RESULTS – 7 / 14 / 28 DAYS

The meter allows you to display average results from different periods:

- 7 days
- 14 days
- 28 days

1. Enter the device memory function.
2. Navigate through the stored information using the device buttons.
3. Select the appropriate period.
4. The screen will display the average calculated based on the saved results.

-
5. The average is used to help observe changes in results over time.
 6. The average function does not replace the individual interpretation of results by a physician

Battery replacement instructions:

- Please insert the batteries correctly, paying attention to the +/- polarity.
- Only one type of battery should be used.
- Always replace the entire set of batteries.
- Exhausted batteries should be removed from the device and never thrown into fire.
- Never try to charge regular batteries.
- Store unused batteries in their packaging away from metal objects.
- Do not leave children unsupervised.

Cleaning and maintenance

- Before cleaning, remove the test strip from the meter.
- Wipe the outer surface of the device with a soft, slightly damp cloth.
- Do not immerse the meter in water or other liquids.
- Do not wash the device under running water.
- Be especially careful not to spill any water on the test strip slot.
- Do not allow blood to enter the device.
- Wipe the display gently without applying strong pressure.
- Keep the lancing device clean.
- Do not leave a used lancet in the lancing device.
- Do not use sharp objects for cleaning.
- Do not use strong solvents or abrasives.
- After cleaning, allow the device to dry completely.
- Store the strips in the designated container.
- The meter should be protected from strong impacts and falls.
- Store your device and accessories in a clean and dry place.

Safety tips and rules

- Before first use, read the instruction manual carefully.
- Use the device only for its intended purpose.
- The glucometer is designed for self-monitoring of blood glucose levels.
- The measurement result does not constitute an independent medical diagnosis .
- Do not change your treatment or medication dosage based solely on a single measurement without appropriate medical advice.
- If you experience any disturbing results, please contact your doctor.
- If the result does not correspond to how you feel, check whether the test was performed correctly and take another measurement.
- Pay attention to the unit of the result – mmol/L and mg/dL are not the same numerical values.
- Wash and dry your hands thoroughly before each measurement.

-
- Do not perform the test with dirty or wet hands.
 - Use only appropriate test strips.
 - Do not reuse test strips.
 - Do not use damaged or dirty strips.
 - Check the expiry date of the strips.
 - Do not leave the strip container open for a long time.
 - Protect your strips from moisture.
 - Lancets are intended for single use only.
 - Use a new lancet for each test.
 - Do not share lancets with other people.
 - Do not touch the lancet tip.
 - Do not point the lancing device at your face or other people.
 - Be especially careful when disposing of the used lancet.
 - Used strips and lancets may have come into contact with blood and should be treated as potentially contaminated material.
 - Used items should be disposed of in a way that minimizes the risk of contact with blood by other people.
 - Do not leave lancets or used strips within the reach of children.
 - Keep the meter, lancing device, strips and lancets out of the reach of small children.
 - Do not allow children to play with the lancing device.
 - Keep the button cell battery out of the reach of children.
 - Swallowing a button cell battery may pose a serious health risk and requires immediate medical attention.
 - Do not put batteries in your mouth.
 - Do not short-circuit the battery terminals.
 - Do not throw batteries into fire.
 - Do not attempt to recharge a battery if it is not designed to be recharged.
 - Install the battery with correct polarity.
 - Protect your device from water.
 - Do not immerse the meter in liquid.
 - Do not allow the belt slot to become flooded.
 - Do not insert metal or sharp objects into the device slot.
 - Do not use a device with a cracked casing.
 - Do not use the meter after a serious fall or if you suspect it is damaged.
 - Do not disassemble or modify the device yourself.
 - Do not store the meter in a place exposed to extreme temperatures or high humidity.
 - Do not leave the device in direct sunlight for long periods of time.
 - Do not place heavy objects on the display.
 - Do not test from your forearm in situations where your glucose level may change rapidly.
 - If the forearm result is questionable, repeat the test with blood taken from your fingertip.

- The memory and average results function is of an auxiliary nature.
- The 7-, 14-, and 28-day averages do not replace medical judgment.
- Do not use the device or accessories if there is any visible damage.
- After each measurement, secure the meter and accessories and store them in a clean, dry and safe place .

Safety instructions for batteries/rechargeable batteries

- Keep batteries/rechargeable batteries out of the reach of children. If swallowed, seek medical attention immediately!
- Disposable batteries must not be recharged. Batteries/rechargeable batteries should not be short-circuited and/or opened. This may result in overheating, fire, or explosion.
- Never throw batteries/rechargeable batteries into fire or water.
- Batteries/rechargeable batteries should never be subjected to mechanical stress.
- Risk of acid leakage from batteries/accumulators.
- Avoid extreme conditions and temperatures that may affect batteries/accumulators, e.g. radiators/direct sunlight.
- If batteries/rechargeable batteries leak, avoid contact with skin, eyes, and mucous membranes! Immediately flush the affected areas with clean water and seek medical attention!
- Leaking or damaged batteries/rechargeable batteries can cause chemical burns when in contact with skin. Therefore, wear appropriate protective gloves.
- The product has a built-in battery that cannot be replaced by the user.

Warning and information symbols

	<u>READ THE INSTRUCTIONS FOR USE</u> Reading the instructions before use will help you perform the measurement correctly, use the meter properly and reduce the risk of incorrect results or incorrect use of the device .
	<u>TIPS AND INFORMATION ON THE MANAGEMENT OF USED PACKAGING</u> The packaging is made of environmentally friendly materials that can be disposed of at your local recycling center. Used packaging material should be delivered to a designated waste disposal site designated by local authorities. For information on recycling options for used products, please contact your local municipality or city office.
	<u>DISPOSAL OF USED ELECTRICAL AND ELECTRONIC EQUIPMENT</u> For environmental reasons, used electrical and electronic products should not be disposed of with regular household waste, but should be disposed of properly. Information on collection points and their opening hours can be obtained from the relevant authorities.

	<u>BIOLOGICAL HAZARD</u> Used test strips, lancets, and other items that come into contact with blood may pose a potential biohazard. Do not reuse them or touch contaminated surfaces. After testing, they should be safely secured and disposed of according to local sharps or biohazard waste management policies.
	<u>CAUTION – SHARP ELEMENTS / STABILITY HAZARD</u> Lancets contain sharp edges that can cause punctures, lacerations, and contact with blood. Exercise extreme caution when inserting, using, and disposing of lancets. Do not touch the blade and do not reuse the lancet. Store the used lancet in a puncture-resistant container and dispose of it according to local regulations.
	This product complies with the requirements of the applicable European and national directives. The product meets European and national requirements for the safety of devices and products.
	This product complies with the requirements of the relevant European and national RoHS directives.
	<u>Improper disposal of batteries/accumulators poses a threat to the environment!</u> Batteries/accumulators should not be disposed of with household waste. They may contain harmful heavy metals and should be treated as special waste. The chemical symbols for heavy metals are: Cd = cadmium, Hg = mercury, Pb = lead. Therefore, used batteries/accumulators should be taken to municipal hazardous waste collection points. The built-in battery must not be disassembled for disposal. The product should be taken to a collection point for used electronics.

We reserve the right to make changes to the text, design and product data without notice.

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